

From Production to QC, Ensure consistency at every drop

Auto titrator

COM-A19 with sample changer



Excellent accuracy

Error free titration with
calculation ready for you

Automatic titration


12 specimen samples
measured in one go

Interchange titrant made easy


Changing the titrant is easy,
just remove the buret head

【Measurement results of total alkalinity in caustic
soda】

No.	Sample Vol (g)	Titrant Vol (mL)	Total alkalinity (%)
1	80.0328	49.589	49.7174
2		49.594	49.7224
3		49.574	49.7023
4		49.582	49.7104
Adv. (%)		49.720	
SD (%)		0.013	
RSD (%)		0.03	

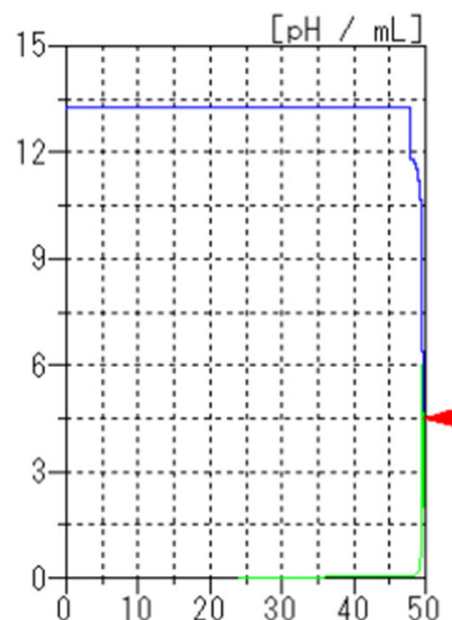

4 easy steps for measurement!

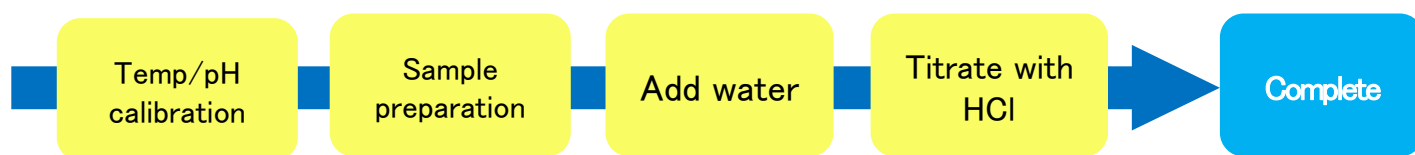


Various configuration for your budget and needs

Affordable



For simple acid/salt analysis

12 specimen sampler

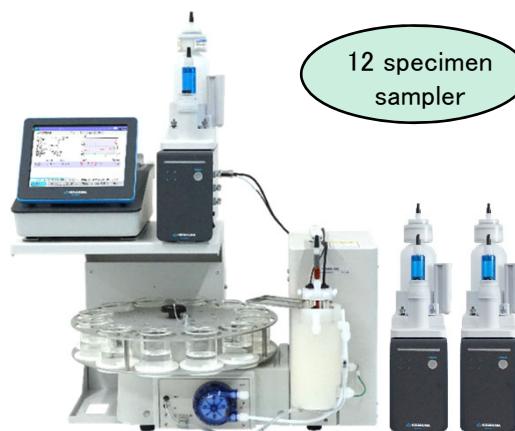

For frequent automatic measurement with lots of samples

8.4" touch panel



High-end model for heavy users, with bigger monitor and built-in thermal printer

24 specimen sampler



Contact us for application data you're looking for!
Application library available from here: